

# ONETOUCH Reveal®



## OneTouch Reveal® web app

Helps your doctor see more about your blood sugar.

- Tracks your blood glucose, food and insulin dosing over time.
- Reports can be printed, ready for discussion at your next appointment with your healthcare professional
- Can be used on PC and Mac computers at [OneTouchReveal.ca](https://www.onetouchreveal.ca)

# ONETOUCH®

# How to setup your OneTouch Reveal® web app

Before you begin, confirm that cookies are enabled in your browser.



Bluetooth

## Wireless Option

OR

## Cable Option

### 1. Connect wireless devices and activate your account

- Install the OneTouch Reveal® mobile app on your mobile device
- Open the mobile app and follow the on-screen instructions to setup your account and pair your OneTouch Verio Reflect® or OneTouch Verio Flex® meter

For more details, refer to the OneTouch Reveal® mobile app **Help File**.

### 2. Sign in to the OneTouch Reveal® web app on your computer

- Go to **OneTouchReveal.ca** and enter the same **Username** and **Password** that you used in the OneTouch Reveal® mobile app
- Click **Sign In** and follow the on-screen instructions.

### 1. Install Data Transfer Tool

- Unplug all OneTouch® USB cables from your computer
- **Sign In** to your OneTouch Reveal® account and click **help** from the top right of any screen
- Click the **Download** button under **Data Transfer Tool** and follow the prompts
- Select the appropriate devices
- Click **Finish**. After installation, you may be prompted to re-start your computer.

### 2. Transfer data from your OneTouch® meter

- Connect your meter to the computer with your OneTouch® USB cable (*cable sold separately*).
- Open the OneTouch Reveal® Data Transfer Tool and select **Detect Now**. When the device is found, a confirmation screen will appear.
- **Sign In** to the OneTouch Reveal® web application
- Click **Confirm** to transfer data. When the transfer completes, click **OK**.




You're all set! You can now access your data with the OneTouch Reveal® web app at [OneTouchReveal.ca](https://OneTouchReveal.ca)

**ONETOUCH®**



## View your blood glucose readings

- Click the **Home** tab to see your blood glucose readings
- You can see your colour-coded blood glucose readings, **High (red)**, **In Range (green)**, **Low (blue)**, over the previous 14 days.

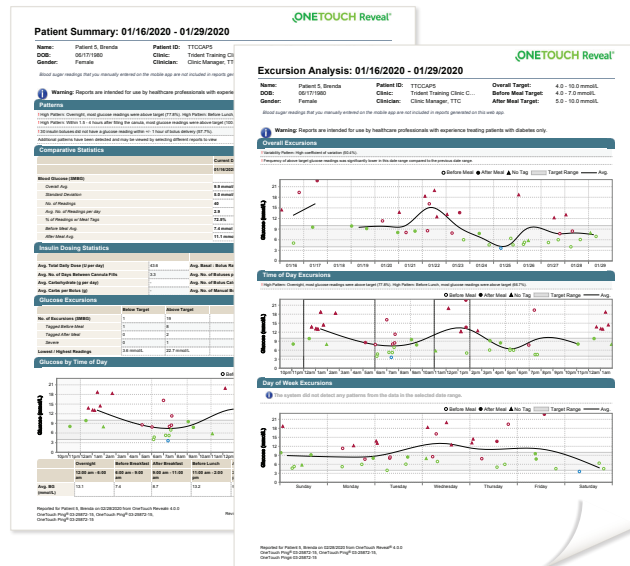

## Prepare for your next visit with your healthcare professional

Select the **Report** you want to discuss with your healthcare professional. Click **Download Report**.

# Any questions? We're here to help



Visit **OneTouch.ca**



Call us at 1-800-663-5521

Find the necessary system requirements at  
[OneTouchReveal.ca/support/en\\_CA/](https://OneTouchReveal.ca/support/en_CA/)

